


BVF™ Oval Filters


Now available with integrated oval mouthpiece

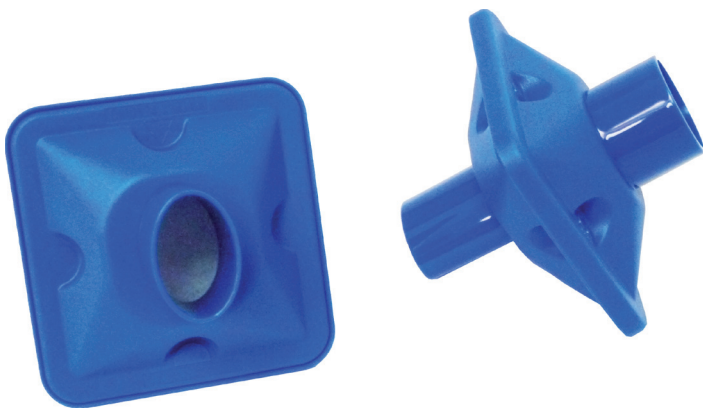

Specially designed in response to the patients need.

Unlike barrier filters which trap expectorated matter whilst allowing viruses and bacteria to pass through, Vitalograph filters use an electrostatically charged material that traps expectorated matter plus bacteria and viruses for highly effective protection against cross contamination.

In addition to the hygiene benefits of using bacterial viral filters they also support productivity in pulmonary function testing by reducing the time spent cleaning and sterilising test equipment.

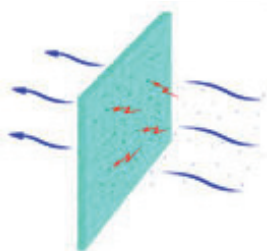

- ★ Single use filters for effective filtration of bacteria, viruses and expectorated matter.
- ★ Suitable for inspiratory and expiratory spirometry testing.
- ★ Protects the subject against cross contamination.
- ★ Protects the spirometer against cross contamination.
- ★ 99.99% cross-contamination efficiency at 1 micron.
- ★ Very low expiratory flow impedance.
- ★ Low effective dead space of 65ml as required for lung volume and DLCO testing.
- ★ Comfortable mouthpiece.
- ★ Efficient at high flow rates.
- ★ Attractively priced.



Pulmonary Function Test Kits

- ★ Complete pulmonary function test kits for use with Jaeger spirometers.
- ★ Kit includes BVF bacterial/viral filter, plastic mouthpiece, elbow and nose clip.
- ★ Completely disposable.



Bite-On Mouthpieces

- ★ Protects the subject from cross infection.
- ★ High quality bite-on mouthpieces designed for use with any size filter.
- ★ Kit includes BVF bacterial/viral filter, plastic mouthpiece, elbow and nose clip.
- ★ Completely disposable.



SafeTway® Mouthpieces

- ★ Protects the subject from cross infection.
- ★ Suitable for expiratory manoeuvres with peak flow meters and spirometers
- ★ One-way valve prevents inhalation respiratory manoeuvre.
- ★ Completely disposable.



Noseclips

- ★ Designed for a comfortable fit.
- ★ Completely disposable.

Colour:	Code with oval mouthpiece	To fit spirometers from:
Blue	28750	Vitalograph, MIR, Micromedical, other office spirometers
White	28762	Sensormedics, Jaeger, Zan, PK Morgan, other Lab PF Equipment
Green	28430	nSpire, ferraris, PDS and KoKo