

Safe Ear Curettes™

Bionix® Safe Ear Curettes allow healthcare providers to safely and gently remove cerumen with styles to match different types of wax and personal preference.



PROBLEM

Hard, rigid stainless steel curettes not only cause patients pain and discomfort, but can also easily damage the ear canal and tympanic membrane.

SOLUTION

Bionix Safe Ear Curettes are flexible, to reduce injury to the ear canal, making them safer than stainless steel curettes.

RESULTS

Seven color coded tip styles to meet specific needs, giving you “the right tool for any job” for every cerumen situation.

*“I have been using Bionix ear curettes since residency.
To me, they are the gold standard.” - Sidney Shah, MD*

POTENTIAL REIMBURSEMENT (U.S. ONLY)

CPT Code/Description

69210 - For potential reimbursement of impacted cerumen.



Obelis s.a • Boulevard Général Wahis 53 • 1030 Brussels, BELGIUM

Tel: +(32) 2.732.59.54 • Fax: +(32) 2.732.60.03 • E-Mail: mail@obelis.net

Bionix LLC • 1670 Indian Wood Circle, Maumee, Ohio 43537, United States

Phone: +1 419.727.8421 • Fax: +1 800.455.5678 • www.Bionix.com



INNOVATIVE SOLUTIONS FOR EVERYDAY PROBLEMS

800.551.7096 • www.bionix.com

©2021 Bionix LLC | MKT-0045 | Rev. D

Safe Ear Curettes™

PRODUCT OFFERING

ADULT

Yellow CeraSpoon® | #6333
The Yellow CeraSpoon is perfect for softer wax, especially in patients who produce excessive amounts of cerumen. **4mm**

Orange ControLoop® | #2999
The Orange ControLoop has a bent “Shapleigh” style tip that allows the wax to be pulled straight out of the ear canal. **4mm**

Red AngleLoop® | #3444
The Red AngleLoop has an angled tip that facilitates reaching the hard-to-get, tough cerumen common in older patients. **4mm**

GENERAL

White FlexLoop® | #9555
The White FlexLoop is an all-purpose curette designed for normal curettage in a variety of patients. **4mm**


YOUTH

Purple VersaLoop® | #4111
The Purple VersaLoop has a smooth, flexible tip that gives the user exceptional control while offering the patient more comfort. **3mm**

Blue InfantScoop® | #4888
The Blue InfantScoop is designed for smaller ears and may be used through an otoscope. **2mm**

Green MicroLoop® | #1222
The Green MicroLoop is more rigid in construction and is designed to remove tough, impacted cerumen. **3mm**

USAGE INSTRUCTIONS

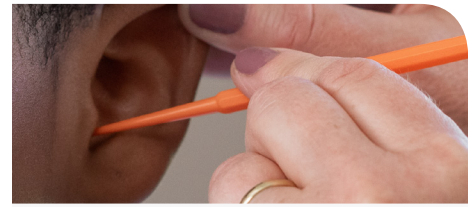

NORMAL CURETTAGE

1. Use otoscope to determine cerumen location.
2. Choose the appropriate curette for the type of patient and cerumen.
3. Place curette into ear canal and remove cerumen.
4. Repeat otoscope exam and determine location of any remaining cerumen.



OTOSCOPE CURETTAGE

1. Use otoscope to determine cerumen location.
2. Select one of the Youth Safe Ear Curettes™.
3. Slide scope magnification to one side.
4. Place curette through otoscope and remove cerumen.

“The variety of sizes and tip shapes makes ear wax removal easier than ever before.” - Carl A. Moritz, Jr., MD

Refer to the package insert supplied with product for specific information on indications for use, contraindications, warnings, and precautions.

BIONIX®, CeraSpoon®, FlexLoop®, InfantScoop®, MicroLoop®, AngleLoop®, VersaLoop® and ControLoop® are registered trademarks. Safe Ear Curettes is a trademark of BIONIX.



INNOVATIVE SOLUTIONS FOR EVERYDAY PROBLEMS

800.551.7096 • www.bionix.com